

L.DROP

Dynamic. Accurate. Connected.

Non-contact Dispensers




L.DROP

A contemporary liquid handler for reproducible genomics and drug discovery.

The L.DROP Liquid Handler provides users with a connected experience that is unrivaled by traditional liquid handling instruments.

The platform is ideal for laboratories looking to leverage state-of-the-art automation to increase productivity, reduce errors, streamline workflows, optimize assay performance, and eliminate repetitive tasks.

- Vision-based deck-calibration for a quick and easy installation.
- Cloud-based protocol generation allows any researcher to program complex methods and collaborate with others regardless of automation experience.
- Automatic liquid class creation enables a greatly simplified optimization and implementation process for protocols.
- Dynamically verified volumetric accuracy and precision for a wide range of samples, reagents, and buffers guarantees the generation of trustworthy data
- Vision-based labware detection ensures that users cannot make mistakes placing labware and setting up runs
- On-demand sample traceability provides visibility to the sample and reagent contained in any well at any point in the protocol.
- Thoughtful error recovery options protect precious samples and costly reagents from waste.

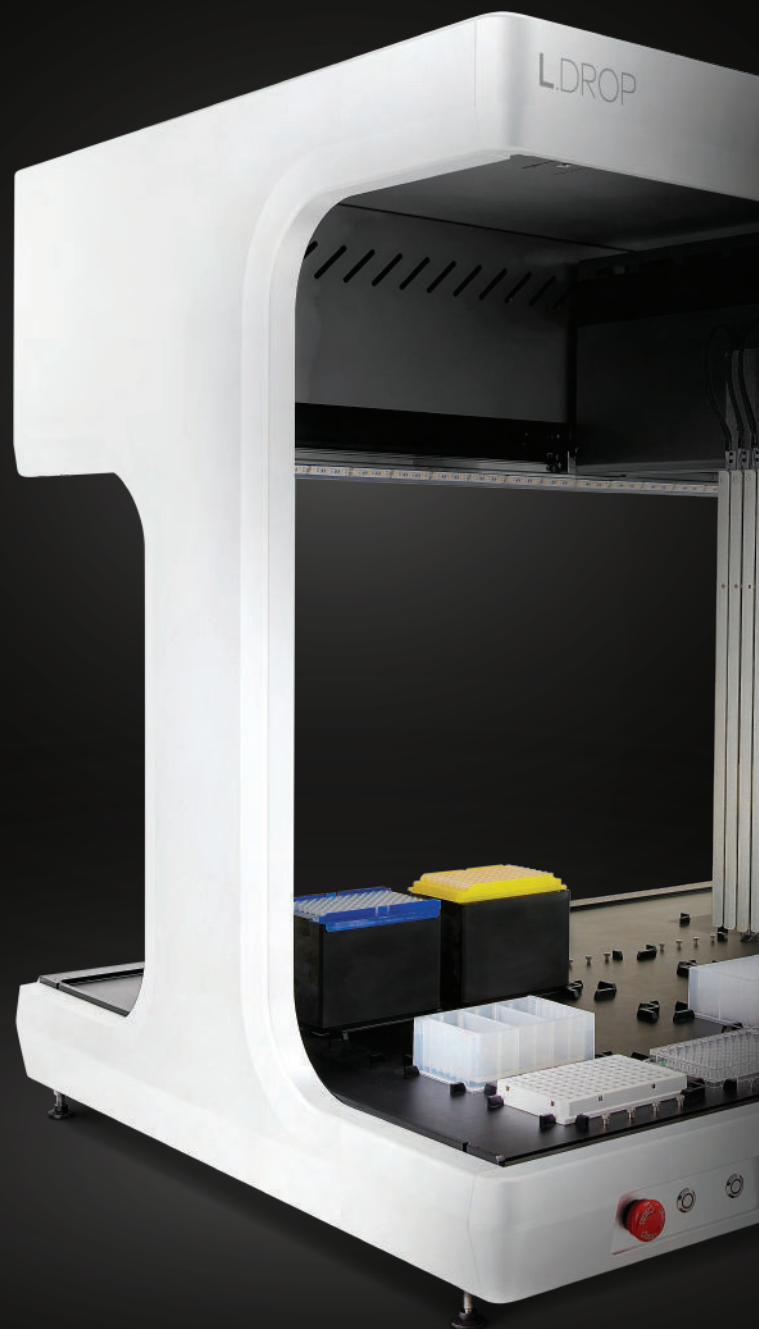

Key Benefits



Swappable 8-channel pipetting head with a volume range from 1 – 1000 μ L.



Pressure and flow Sensors for Real-time Volume Verification.



Zero-force tip attachment for reliable and wear-free hardware.



Camera with Cloud-based Image Analysis for Deck Verification.



High-resolution Camera for Automated Optical Calibration.



Real-time deck monitoring to capture video of errors for future troubleshooting and correction.



Algorithms for Automated Liquid Class Creation.



Holds up to 20 SBS format microplates, tips, or lids for easy setup and switching between protocols and five integration positions for on-deck thermal control and shaking modules.



Optimized to handle diverse applications.

DISPENDIX has developed a growing library of protocols which are accessible on-demand from the cloud platform.



CELL AND GENE THERAPY /SYNTHETIC BIOLOGY

- Cell culture
- CRISPR Gene Editing
- Protein Characterization
- Protein and Nucleic Acid Quantification



GENOMICS

- Library Preparation for DNA-Seq, RNA-Seq, and Targeted Sequencing
- qPCR Plate Preparation



DRUG DISCOVERY

- Assay Development and Optimization (DoE)
- Immunoassays
- Cell Based Assays
- Biochemical Assays
- High-throughput Screening



GENERAL LIQUID HANDLING

- Serial Dilution
- Reagent Formulation and Aliquoting
- Sample Reformatting
- Cherry Picking
- Normalization and Pooling



Productivity from the cloud down.

The L.DROP is run through an intuitive Cloud-based platform that enables users to build, run, and analyze protocols from anywhere – anytime!

1

Build It

- Build protocols using a straightforward drag and drop interface.
- Program methods that are focused on the scientific method rather than the automation.
- Automatically "fit" protocols to liquid handlers and laboratories of varying configurations and resources.
- Capture instrument specific constraints and iterate rapidly with on-demand simulations rather than trial and error wet-runs.

2

Run It

- Run and optimize protocols from the cloud with the first ever liquid handling experience that does not require you to be tethered to the instrument.
- Monitor workflow status, instrument status, errors, and protocol execution from the comfort of your office or even from home.

3

Analyze It

- Access data for any sample retroactively through state and liquid tracking information from the cloud.
- Optimize laboratory logistics with visibility to, reagent and consumable consumption.

Connected automation brings productivity and scale to your lab



Spanning Channels
for Transfers Across
Varying Labware

Air Driven Pipetting
from 1-1000 μ L with
Integrated Flow and
Pressure Sensors

Power and
Emergency Stop

Zero-force Pipette
Tip Attachment

Through-deck
Tip Chute

Seamless Enclosure
Fits in Standard
Hood


Swappable
8-Channel
Pipetting Module

Built-in Positions
for Thermal and
Shaking Accessories

Finally, a user experience focused instrument for liquid handling automation

L.DROP offers cutting-edge technology in automation to enable researchers to program and execute automated liquid handling methods regardless of automation experience.

- Automate deck calibration and liquid classes with support from intelligent sensors
- Build protocols and collaborate with your colleagues from anywhere – anytime!
- Run and aggregate data from all the instruments in your workflow from a single cloud-based platform.
- Generate trustworthy data with real-time volume verification and sample traceability in the cloud.

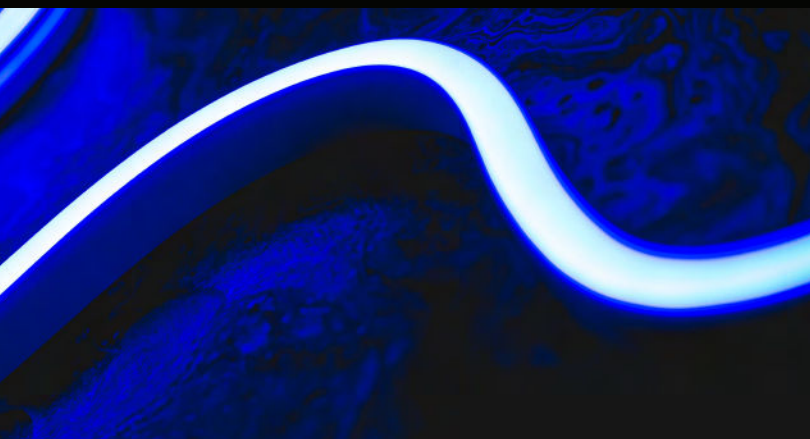

L.DROP



©2021 BICO AB. All rights reserved. Duplication and/or reproduction of all or any portion of this document without the express written consent of BICO is strictly forbidden. Nothing contained herein shall constitute any warranty, express or implied, as to the performance of any products described herein. Any and all warranties applicable to any products are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. BICO provides no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein. The use of products described herein is subject to certain restrictions as set forth in the applicable terms and conditions of sale accompanying the purchase of such product. BICO may refer to the products or services offered by other companies by their brand name or company name solely for clarity and does not claim any rights to those third-party marks or names. BICO products may be covered by one or more patents. The use of products described herein is subject to BICO's terms and conditions of sale and such other terms that have been agreed to in writing between BICO and user. All products and services described herein are intended FOR RESEARCH USE ONLY and NOT FOR USE IN DIAGNOSTIC PROCEDURES.

The use of BICO products in practicing the methods set forth herein has not been validated by BICO, and such nonvalidated use is NOT COVERED BY BICO'S STANDARD WARRANTY, AND BICO HEREBY DISCLAIMS ANY AND ALL WARRANTIES FOR SUCH USE. Nothing in this document should be construed as altering, waiving or amending in any manner BICO's terms and conditions of sale for the instruments, consumables or software mentioned, including without limitation such terms and conditions relating to certain use restrictions, limited license, warranty and limitation of liability, and nothing in this document shall be deemed to be Documentation, as that term is set forth in such terms and conditions of sale. Nothing in this document shall be construed as any representation by BICO that it currently or will at any time in the future offer or in any way support any application set forth herein.

